

Clario



MOBILITY AND RELIABILITY



Clario Airway Suction Pump

Clario

Easy to use and highly mobile suction

Home healthcare teams and patients in need of airway suctioning appreciate the Clario suction pump because of its clever engineering and smart design.

Clario's modular design consists of 3 main parts for easy, uncomplicated assembly:

- Drive unit
- Canister
- SafetyChamber

Clario has 3 preset vacuum levels that make handling easy.

Clario is very hygienic as no secretions can enter the pump, making routine cleaning simple.

Clario is quickly ready for the next user.

Clario's lightweight design means it can be taken everywhere. Clario is available in a battery-operated version for even more mobility.



- a** Floater prevents an overflow
- b** Plugs for safe transport
- c** 0.5 L canister
- d** Vacuum regulator
- e** SafetyChamber
- f** LED for battery status
- g** Carrying handle

High mobility combined with high reliability

Clario combines clever engineering with smart design and ease of use. It consists of 3 components only, all of which can be clicked together in an instant. The tubing required for suctioning is conveniently attached to the canister lid. Clario offers 3 easy-to-set, easy-to-explain vacuum settings. An LED display on the AC/DC version indicates when battery should be charged. All these features are presented in an appealingly non-medical, discreet design and at a surprisingly low noise level.



Easy to use

Threefold hygienic design

Clario has a threefold hygienic design. The floater prevents an overflow. The SafetyChamber will collect secretions if the floater was forgotten. Finally, the membrane prevents any secretions from entering the pump. These measures ensure that there is no overflow into the mechanical parts of the pump. There is no need to open the device for cleaning and no filters are required. The pump is very easy to clean to help prevent cross-contamination amongst patients.



Easy to clean

Compact and portable

Clario offers more freedom for the user. The battery-operated version, weighing only 2 kg, features a rechargeable battery that allows for more than 50 minutes of suction time. The optional carrying bag for the Clario and its components makes it discreet and convenient to transport.



Easy to transport

Fast Facts

- Medela offers healthcare professionals 50 years of experience in medical vacuum technology.
- Our swiss development and manufacturing process yields mechanical precision to guarantee excellent performance.
- Clario is appreciated in the homecare setting for reliability and ease of use.
- Clario offers modern yet timeless design. It does not look like a typical medical device.
- Clario is highly recommended by Home Healthcare Teams.
- Clario needs little time and effort for cleaning.
- Clario's threefold hygienic design helps eliminate cross-contamination.
- Clario increases mobility, thus improving the quality of life.
- 2 years warranty on device.

Accessories

Plug for 12 VDC IN car connection cable
for operation in a car.



Vacuum gauge
for controlling vacuum settings.


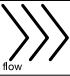



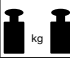
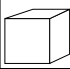
Soft case carrying bag
for simple and discreet transport.






Use Medela accessories for even greater flexibility and ease of use.

Technical data

| | | |
|-------------------------------------------------------------------------------------|----------------------------------------------|-------------------------------------|
|  | min. -18 kPa med. -36 kPa max. -80 kPa | -135 mmHg -270 mmHg -600 mmHg |
|  | low flow 15 L/min. | |
|  | Power 100 - 240 V 50 - 60 Hz | |

| | | |
|-------------------------------------------------------------------------------------|--------------------------------|------------------|
|  | AC AC/DC | 1,7 kg 2,0 kg |
|  | H x W x D 223 x 255 x 95 mm | |

| | | |
|---------------------------------------------------------------------------------------|------------------------------------------------------------------|---------------------------------------------------------------------------------------|
|  | ISO 9001 ISO 13485 CE (93/42/EEC), IIa 2 years warranty |  |
|  | | |

Medical Vacuum Technology for Healthcare Professionals

Please contact us or your local Medela representative for details.

 Medela AG
Lättichstrasse 4b
6341 Baar, Switzerland
www.medela.com

Local contact:

Medela Australia Pty Ltd
Medical Technology
3 Arco Lane
Heatherton, Vic, 3202
Australia
Phone +61 (0) 3 9552 8600
Fax +61 (0) 3 9552 8699
contact@medela.com.au
www.medela.com.au

Medela Canada Inc.
4160 Sladeview Cres., Unit #8
Mississauga, Ontario
L5L 0A1
Canada
Toll Free 800.435.8316
Toll Free Fax 800.995.7867
info@medela.ca
www.medela.ca

Medela UK
Huntsman Drive
Northbank Industrial Park
Irlam, Manchester M44 5EG
UK
Phone +44 870 950 5994
Fax +44 870 389 2233
info@medela.co.uk
www.medela.co.uk

Medela Inc. USA
1101 Corporate Drive
McHenry, IL 60050
USA
Phone +1 877 735 1626
Fax +1 815 363 2487
suction@medela.com
www.medelasuction.com