


OCULUS Maddox-Kreuz/ Maddox cross

Tangentenskala/Tangent scale

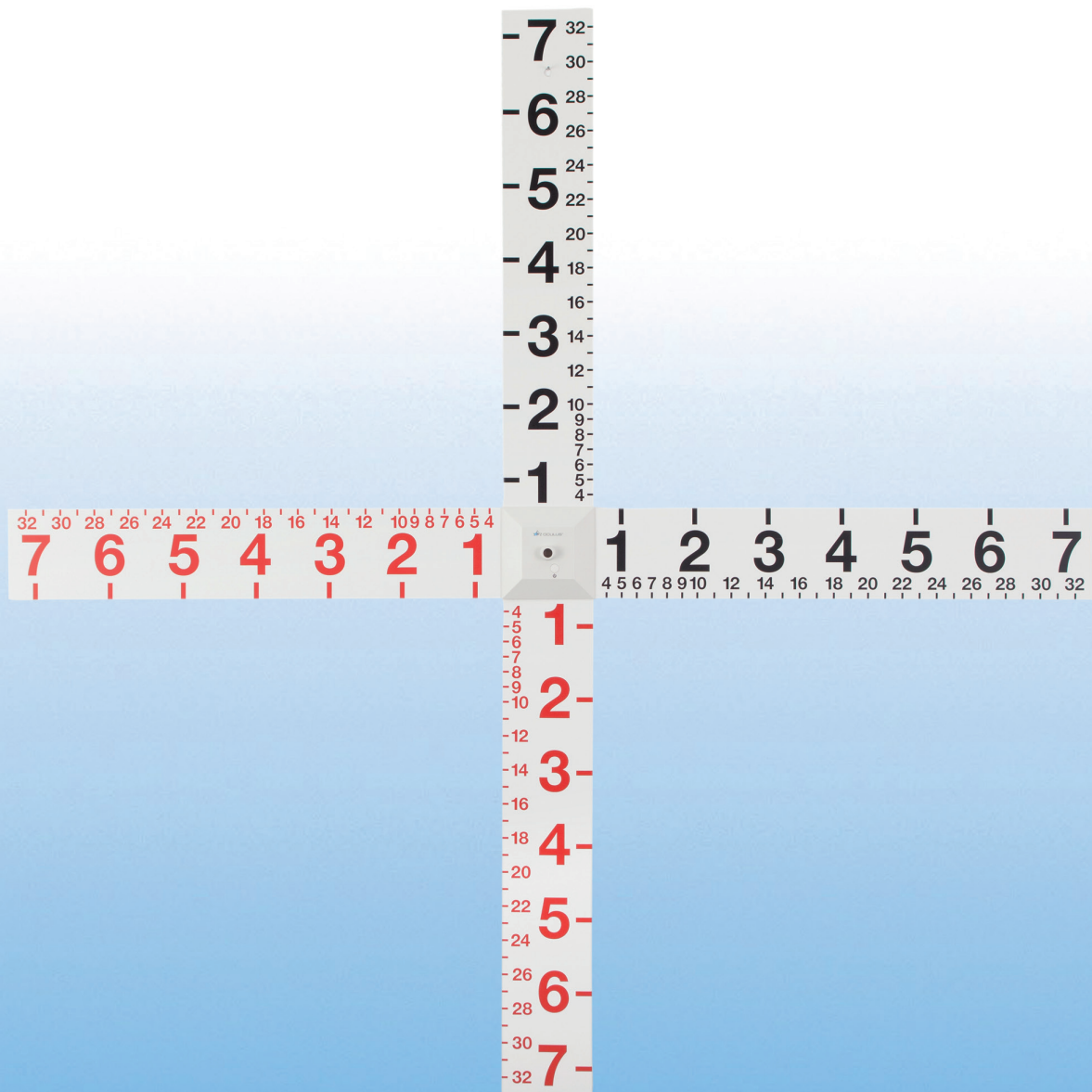

OCULUS Maddox-Kreuz/ OCULUS Maddox cross

Tangentenskala zur Schielwinkelbestimmung und Korrespondenzprüfung

Das Maddox-Kreuz besteht aus einem Skalenset mit zwei unterschiedlichen Skalierungen:

- große Ziffern für 5 m Prüferentfernung
- kleine Ziffern für 1 m Prüferentfernung

Die Zahlen geben die Winkelgrade an. Im Zentrum der Skala befindet sich ein punktförmiges LED-Fixierlicht, das sich bei der Untersuchung in Augenhöhe genau gegenüber dem Patienten befinden sollte.



> Lieferumfang/Delivery contents

Tangent scale for squint angle measurement and retinal correspondence testing

The Maddox cross is provided with two scales:

- large-digit scale for 5 m testing distance
- small-digit scale for 1 m testing distance

The numbers denote the angle in degrees. In the centre of the scale is a spot fixation light which should be located at eye level exactly opposite the patient during the examination.

Technische Daten/Technical Data

Technische Angaben		Technical specifications	
Abmessung (B x T x H)	1290 x 28 x 1290 mm	Dimensions (W x D x H)	58.8 x 1.1 x 58.8 in
Gewicht		Weight	
Maddox-Kreuz	920 g	Maddox cross	2.0 lbs
Netzadapter	65 g	Power adapter	0.15 lbs
Stromversorgung	5 V DC, 100mA, USB-Stecker Typ A	Power supply	5 V DC, 100mA, USB plug type A
Mitgeliefertes Steckernetzteil	USB, 100-240 V AC, 50-60 Hz 0,3 A, Sek. 5 V DC, 1,2 A, REF 05150950, inkl. EU-Adapter	Provided power adapter	USB, 100-240 V AC, 50-60 Hz 0,3 A, sec. 5 V DC, 1,2 A, REF 05150950, incl. EU adapter
Beleuchtung	LED	Lighting	LED

CE gemäß Richtlinie 93/42/EWG über Medizinprodukte/in accordance with Medical Device Directive 93/42/EEC

WWW.OCULUS.DE



OCULUS ist zertifiziert gemäß
DIN EN ISO 13485/
OCULUS is certified by TÜV according to
DIN EN ISO 13485

OCULUS Optikgeräte GmbH
Postfach • 35549 Wetzlar • GERMANY
Tel. +49-641-2005-0 • Fax +49-641-2005-295
E-Mail: sales@oculus.de • www.oculus.de

- OCULUS USA, info@oculususa.com
- OCULUS Asia, info@oculus.hk
- OCULUS Czechia, oculus@oculus.cz
- OCULUS Iberia, info@oculus.es
- OCULUS Poland, biuro@oculus.pl
- OCULUS Turkey, info@oculus-turkey.com.tr