

OCULUS Ophthalmic Tables

Universal lifting tables
for medical instruments



OCULUS Ophthalmic Tables

High-quality and ergonomic

OCULUS Ophthalmic Tables catch the eye with their distinctive design. Their multiplex tabletops with state-of-the-art high pressure laminate (HPL) give them a classy, aesthetic appearance. The HPL coating is highly resistant to chemical and thermal influences. The smooth surface is characterized by extreme strength and durability – making it just perfect for use in the medical field.



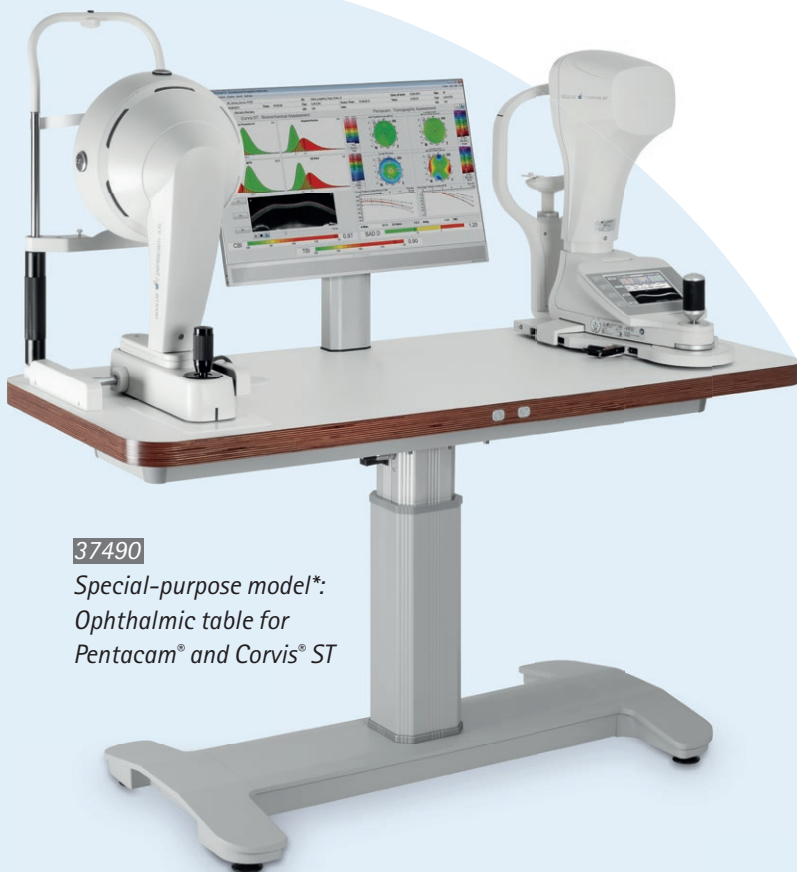
37420 *Ophthalmic table for one device, symmetrical*



37430 *Ophthalmic table for one device, asymmetrical*



37440 *Ophthalmic table for one device, asymmetrical, wheelchair-accessible*



37490
Special-purpose model:
Ophthalmic table for
Pentacam® and Corvis® ST*



37450 *Ophthalmic table for two devices, delta-shaped tabletop*



37460 *Ophthalmic table for two devices, symmetrical*

37475 *Identical to 37460, but with*

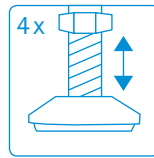
37485 *different tabletop sizes*

* Further special-purpose tables are available for various OCULUS devices. Please contact your local distributor for details.


Designed to fit anywhere in the room

All four feet are height-adjustable, providing a safe and steady base also on uneven floors. The table base is optionally also available on caster wheels.

All ophthalmic tables are equipped with a smooth-sliding vertically adjustable column.



height-adjustable

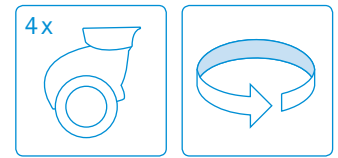

caster wheels with fixing brakes



Illuminated height adjustment keys

The height adjustment keys are smoothly incorporated into the table edge. They are illuminated and are thus easy to find and operate also in darkened rooms.



Cables neatly stowed away

Running from two connectors in the table base up to the tabletop are two network cables for your network devices. Neatly hidden away: The network cables run through the column.



The power connections for your devices are located under the tabletop. Two device connection cables are included.

Tabletop also available in custom colour design



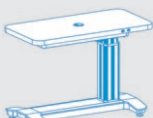

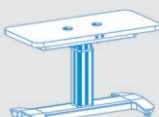

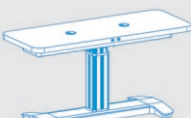
The tabletop comes in discreet matte white. The edges are stained in dark brown. Custom colours on request.

The table foot is powder-coated in light grey.



OCULUS Ophthalmic Tables

Technical Data

	Art. no.	Tabletop	Lift	Weight	Max. load capacity	Electr. connection*
	37420	630 x 500 mm (24.8 x 19.7 in)	655 – 905 mm (26.4 – 34.3 in)	29 kg (63.9 lbs)	84 kg (185.2 lbs)	230 V
	37430	770 x 500 mm (30.3 x 19.7 in)	655 – 905 mm (26.4 – 34.3 in)	33 kg (72.8 lbs)	84 kg (185.2 lbs)	230 V
	37440	980 x 500 mm (38.6 x 19.7 in)	655 – 905 mm (26.4 – 34.3 in)	39 kg (86.0 lbs)	84 kg (185.2 lbs)	230 V
	37450	1050 x 570 mm (41.3 x 22.4 in)	655 – 905 mm (26.4 – 34.3 in)	36 kg (79.4 lbs)	84 kg (185.2 lbs)	230 V
	37460	1050 x 500 mm (41.3 x 19.7 in)	655 – 905 mm (26.4 – 34.3 in)	42 kg (92.6 lbs)	84 kg (185.2 lbs)	230 V
	37475	770 x 500 mm (30.3 x 19.7 in)	655 – 905 mm (26.4 – 34.3 in)	33 kg (72.8 lbs)	84 kg (185.2 lbs)	230 V
	37485	1250 x 500 mm (49.2 x 19.7 in)	655 – 905 mm (26.4 – 34.3 in)	47 kg (103.6 lbs)	84 kg (185.2 lbs)	230 V

* All tables are also available with 115 V.

CE in accordance with Medical Device Directive 93/42/EEC

WWW.OCULUS.DE



OCULUS is certified by TÜV according to
DIN EN ISO 13485 MDSAP

OCULUS Optikgeräte GmbH
Postfach • 35549 Wetzlar • GERMANY
Tel. +49 641 2005-0 • Fax +49 641 2005-295
Email: export@oculus.de • www.oculus.de

Find your local OCULUS representative on our website.