

HEADQUARTERS

Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5568
F +886 2 2268 6525
marketing@apexmedicalcorp.com

PAN-ASIA

TAIWAN
Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5600 /0800 054 133
F +886 2 2268 6526
customer@apexmedicalcorp.com

STURDY INDUSTRIAL

168, Sec.1, Zhongxing Road, Wugu District,
New Taipei City 24872, Taiwan
T +886 2 2981 5245
F +886 2 2984 8208
foreign@sturdy.com.tw
www.sturdy.com.tw

CHINA

Apex Medical (Kunshan) Corp.
No. 1368, ZiZhu Road, Kunshan,
Jiang Su Province, China
T +86 512 8617 1660 Ext 1150
F +86 512 8617 8663
customer.cn@apexmedicalcorp.com

INDIA

Apex Medicalcorp India Pvt. Ltd.
303, Third Floor Block 4A DLF
Corporate Park Phase-III,
Gurgaon, Haryana
T +91 11 6654 4138
F +91 11 6654 4052
enquiries.india@apexmedicalcorp.com

THAILAND

Apex Medical Thailand
65/100 Chamnanphenjati Business Center Bldg.
11/F Rama 9 Road, Huaykwang
Bangkok 10320
T +662 948 4848
F +662 948 4170
sales.thai@apexmedicalcorp.com

MIDDLE EAST • AFRICA

29, Misr Liltameer Building,
Stree#10, AREA#8, Sheraton Building
Heliopolis Cairo, Egypt.
T +202 2 2671303
F +202 2 2678649
raafat.wahby@apexmedicalcorp.com

PAN-EUROPE

SPAIN
Apex Medical S.L.
Elcano 9, 6ª planta 48008
Bilbao, Vizcaya, Spain
T +34 94 4706408
F +34 94 4706409
info@apexmedical.es

UK

Apex Medical Limited
Unit 33, Great Western Business Park
McKenzie Way, Worcester, WR4 9GN
United Kingdom
T +44 (0) 190 5774 695
sales@apexmedicalcorp.co.uk

FRANCE

Apex Medical France
4 Boulevard de la Chanterie
49124 Saint-Barthelemy d'Anjou
France
T +33 2 72 79 27 13
F +33 2 85 73 07 67
info.france@apexmedicalcorp.com


GERMANY

SLK Vertriebsgesellschaft GmbH
Lindenhorster Straße 38-40 44147
Dortmund
T +49 (0)231-925360-175
F +49 (0)231-925360-29
info@slk-gmbh.de
www.slk-gmbh.de

AMERICAS

USA

Apex Medical USA
927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
sales.usa@apexmedicalcorp.com

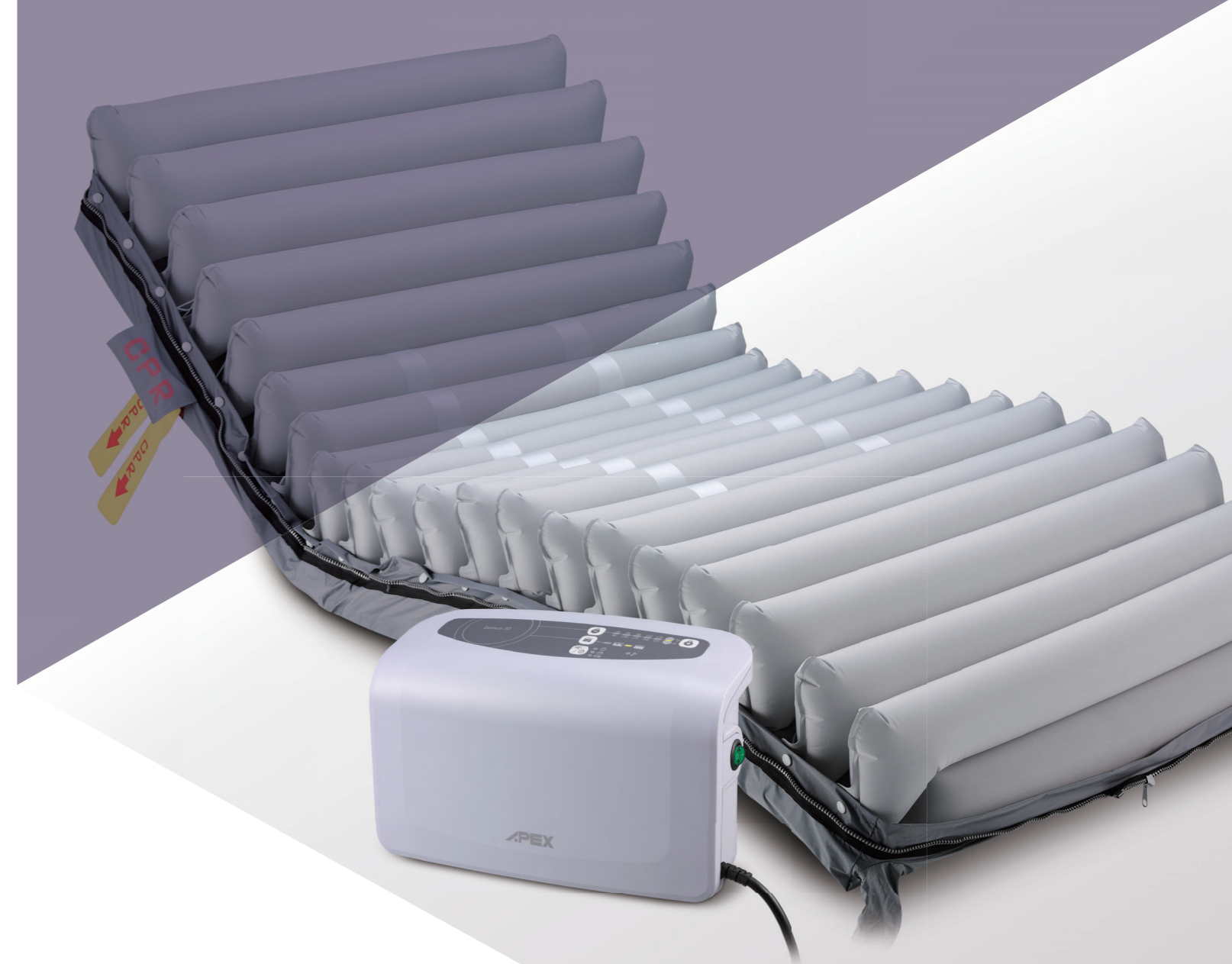

Domus 3D

Alternating Pressure Redistribution System

- Suitable for patients in medium to high risk of pressure injury
- Digital pump for precise pressure control
- Key functions to satisfy pressure ulcer therapy and overall prevention
- Continuous Low Pressure Therapy, increased comfort for better care




Distributed by

Domus 3D


Alternating Pressure Redistribution System

"Your essential solution"

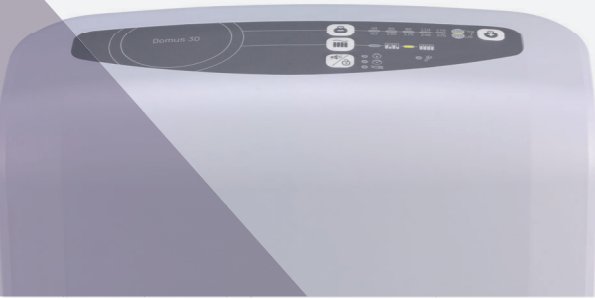

Pressure Relief During Alternating

Domus 3D is an essential solution which not only offers alternating pressure redistribution, but also has a static mode with continuous low pressure

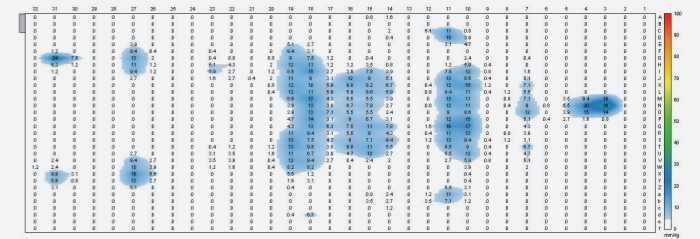

Cell-on-cell design (8")
The bottom cells are capable of providing adequate support for the system without an extra foam base in the case of power failure



Low Air Loss
Provides good ventilation and reduces the accumulation of heat and moisture



Quick Connector with Cap
Provides a stable surface for patient when transporting




CPR Strap
Quickly deflates the mattress in case of an emergency. Designed to prevent accidental trigger



Facile foam compartment (5")
Extra foam interlayer with buttoned design as an adequate support base



Storage strap
Side buckle can simply wrap and store mattress

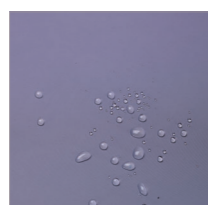

Nylon/PU Cover
Vapour permeable cover offers you the comfort and protection you need



Cable management
Storage strap on the side of mattress for managing cable

Visual and audible alarm


- Power failure alarm
- Low pressure alarm
- Tech. Support alarm

Panel lock
Helps prevent unauthorized changes to the pump settings



2 Modes to Fulfill Individual Needs

- Alternating
- Constant low pressure (CLP)

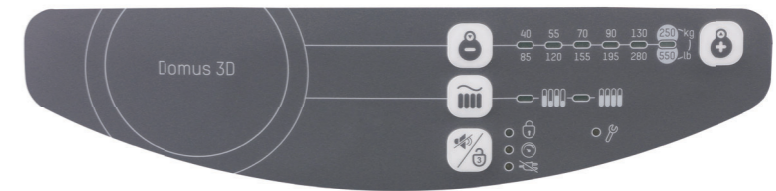

Digital for precision

- Digital pump makes precise pressure control possible for better patient care. Enables the compressor to work more efficiently, creating a longer-lasting system

Alternating

- Continuously and sequentially inflate and deflate air cells to mimic natural sleep movement patterns and avoid long term pressurization of tissue, while stimulating blood circulation

Constant low pressure mode (CLP)

- Reduces the pressure precisely by only 1/3 from the alternate pressure level, providing a static comfort surface while still retaining sufficient support

Specifications

Model	Domus 3D		
Pump	Dimension	29x18.5x12.6 cm	
	Weight	2.3Kg	
	Case Material	Flame retardant ABS	
Supply Voltage	AC 220 - 240V / 50Hz ; AC 120V / 60Hz		
	Mattress	Overlay	Replacement
		Dimension	190 x 85 x 4"
Cell Height	200 x 80/85/90 x 5"+3		
	Weight	17 x 4"/5" cells	20 x 8" cells (cell-on-cell)
	Weight	5.6 kg	6.9 kg
Cover Material	Nylon PU (Quilted)	Nylon PU (Quilted)	
	4 way Stretch (optional)	4 way Stretch (optional)	
Cell Material	Nylon PU	Nylon PU	
Maximum Patient Weight	200 kg	250 kg	
Flame Retardant Standards	EN597-1; EN597-2	EN597-1; EN597-2	

Specifications are subject to change for improvement without notice

